



CHIRAL APPLICATION NOTE

Sample name: AQC-Tyrosine

Chromatogram:

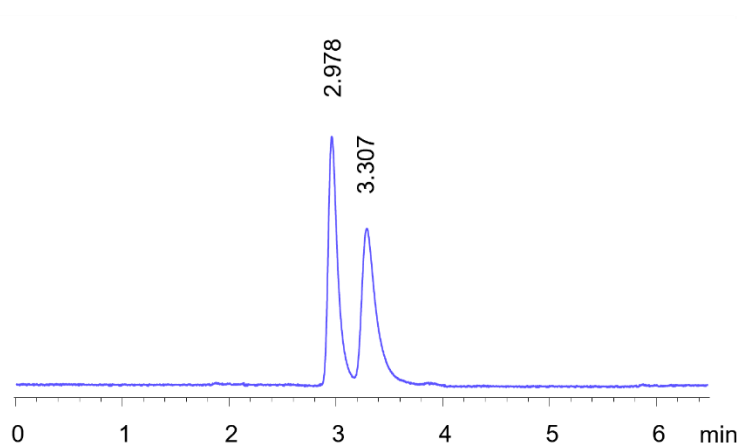

Structure:

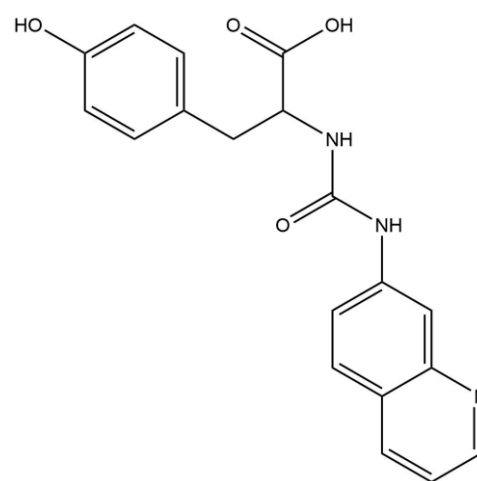

Run Conditions:

Column: TeicoShell (10cmx4.6mm, 2.7 μ m)
Mobile phase: 50/50 Acetonitrile/20 mM Ammonium Formate pH 4.0
Flow Rate: 0.5 mL/min
Temperature: Ambient (23 $^{\circ}$ C)
Injection Volume: 5 μ L
Detection: UV 248 nm
Retention (min): 2.978/3.307 (Peak1/Peak 2)

Resolution: 1.70